Table of contents

1. Warnings	2
2. Contraindications	4
3. Safety features	6
4. Get to know your device	7
4.1 Device description	7
4.2 Intended use	7
4.3 Package contents.....	8
4.4 Device parts.....	8
4.5 Device controls & settings.....	9
4.6 Charging the battery	10
5. Treatment	10
5.1 Treatment Areas.....	10
5.2 Treatments	11
5.3 What to expect	12
5.4 Possible side effects.....	12
6. Get started!	13
7. Cleaning, maintenance & storage	14
8. Troubleshooting	15
9. Customer service	16
10. Warranty	16
11. Disposal	16
12. Signs & symbols	17
13. Specifications	18
14. Legal notice	19

1. Warnings

This device is powerful and should be used with special attention to safety. Before first use, please read this user manual in full, and keep it for future use. Be sure to read all warnings, contraindications and safety precautions, and follow them strictly. We recommend you refamiliarize yourself with the instructions before each treatment.

The device must not be used by the following people: children, transgender women, adults with limited physical, sensory or psychological capacities, or those lacking experience or knowledge in how to use the device in a safe way, or those who do not understand the hazards involved.

- Always ensure that your hands are dry before operating the device, touching the plug, or connecting the charging cradle to the mains.
- Disconnect the device from the mains after charging the battery.
- Keep the device, all parts and accessories away from heated surfaces.

- Never leave the device unattended when plugged into the mains.
- Do not attempt to modify, open or repair the device. You may expose yourself to dangerous electrical components and may cause serious harm to yourself. It will also void your warranty.
- Do not use if the parts or accessories are not supplied by the manufacturer as defined in this manual, or if they are damaged in any way, do not appear to be working properly, or if you see or smell smoke. In such cases, stop using the device, and contact customer service.
- Do not clean any of the parts in the dishwasher.
- Children shall not play with the device. Cleaning and user maintenance shall not be made by children.
- Children must not play with the packaging. They may injure themselves with the packaging material, or swallow small parts and choke.
- This device contains batteries that are non-replaceable.
- This unit has been designed and tested to meet the requirements of electromagnetic, electrostatic, and radio frequency inter-

ference standards. However, the possibility of electromagnetic or other interference may still exist. Relocating the device may help to eliminate interference.

- Do not dispose of the device or the battery with household waste! The law requires that you dispose of this electrical device (including all accessories and parts) at a designated recycling collection point for electrical devices.
- Dispose of the battery safely.
- If the product reaches its end of life, the battery must be removed before disposal (see section 11. Disposal).

2. Contraindications

Do not use if you:

- have a pacemaker or internal defibrillator, or any other active electrical implant anywhere in the body.
- suffer from poorly controlled endocrine disorders, such as diabetes and thyroid conditions.
- have a history of skin cancer, or current condition of any other type of cancer, or pre-malignant moles.
- have received radiation therapy or chemotherapy treatments within the past 3 months, or have areas of potential skin malignancies.

- have an impaired immune system due to immunosuppressive diseases such as AIDS and HIV, or use immunosuppressive medications.
- have used Isotretinoin (Accutane®) within 6 months prior to treatment.
- have a disease related to photosensitivity, such as porphyria, polymorphic light eruption, solar urticaria, lupus, etc.
- suffer from epilepsy.
- are pregnant or nursing.

**Do not use on areas that have:**

- a permanent implant such as metal plates and screws, silicone implants, or areas with injected fills or toxins.
- any active condition in the treatment area, such as sores, psoriasis, eczema, rashes, open lacerations, abrasions, or lesions.
- traumatized or injured skin.
- had recent exposure to the sun without protection (wait 1 day); for tanned or sunburnt skin, wait 3 days.
- had any surgery within 3 months prior to treatment.
- tattoos:
- long or thick hair. Use of the device may result in light bruising.

Consult your physician before continue to use, if you have side effects, such as itchy skin within 48 hours of first use.

Consult your physician if you have:

- couperose or rosacea.
- a history of diseases stimulated by heat, such as recurrent herpes simplex in the treatment area.
- If so, you may treat only following a prophylactic regime.
- used physician-prescribed drugs on your skin within the past 6 months.
- a history of skin disorders, keloids, abnormal wound healing, as well as very dry and fragile skin.
- any other condition which, in the opinion of your physician or dermatologist, would make it unsafe for you to be treated.

3. Safety features

The device protects your skin from overheating.

- **Thermal sensor and heat stabilizer:** These measure the skin temperature, stabilize the heat energy and prevent overheating. When the skin temperature reaches 43°C, the stabilizer will automatically stop emitting heat, until the skin temperature lowers enough for the stabilizer to resume heating.
- **Movement sensor:** The sensor will only emit RF heat when you are moving the device across your skin.

4. Get to know your device

4.1 Device description

The device is an anti-aging device for eyes and face which utilizes RF technology.

This technology affects the fine parts of eyes and face, reaching deep into the dermis for collagen restoration, deep tissue tightening, and reduction for eye wrinkles, facial fine lines and nasolabial lines.

The device has a contact surface (#1) with three electrodes (#2) and a control button (#5) with five surrounding indicator lights.

4.2 Intended use

The device can be used not only for anti-aging of eyes, to reduce and alleviate the fine and dry lines, under-eye bags and dark circles, but also on the face to reduce fine wrinkles on the face (such as nasolabia and the corner of mouth).

The device may be used for eye rejuvenation.

The device is intended for self-treatment in the home and is not suitable for commercial purposes.

Do not use for any purpose other than as described in this manual.

4.3 Package contents

- Silk'n FaceTite Mini
- Type-C adapter
- User manual
- Product warranty

4.4 Device parts

Refer to the diagram on page (i) at the front of this manual to see the location of each feature.

- Treatment surface (#1)
- Electrodes (#2)
- LED light (#3)
- 4 energy level indicator lights (#4)
- Control button (#5)
- Battery indicator light (#6)
- Charging port (#7)

4.5 Device controls & settings

Button	Action	Function / Mode
Control button (#5) 	Long press x1 Short press x1 Short press x2 Short press x3 Short press x4 Long press	Starts the device Energy level 1 - light 1 is on Energy level 2 - light 2 is on Energy level 3 - light 3 is on Energy level 4 - light 4 is on Back to level 1 Shuts down the device
4 energy level indicator lights (#4) 	Shows the energy level selected.	There is one light for each energy level. Level 1 is the lowest energy level, and level 4 is the highest.
Battery indicator light (#6) 	Green - blinks Green - constant Orange - blinks	Device charging Device charged Low charge in battery

4.6 Charging the battery

1. Charge the device before first use.
2. Place the Type-C adapter cable into the device, then connect the adapter to the mains electricity supply. The battery indicator light (#6) will blink green. Charge the device for 3 hours before use. You can use a fully charged device for 20-30 minutes.

Hours

First Charge


Minutes

Regular use

- i** The duration of device use may vary according to different energy levels. Whether the battery is fully charged before use will also affect the service time.

5. Treatment

5.1 Treatment Areas

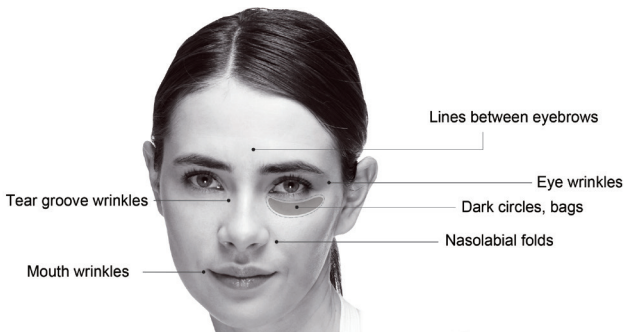5.2 Treatments

- Start with the energy level 1, and change to other energy level after feeling adapt to the energy on the treatment area.
- Use the device for 3 to 5 minutes on each subtle area for maximum effectiveness.
- Use the device 2 to 3 times a week for the first eight weeks.

After eight weeks, use the device weekly, to maintain the effect.

For maximum effectiveness:

- Make sure you apply sufficient treatment medium on the treatment area before use.
- Do not keep the device stand still on one area. Keep moving the device covering the entire treatment area. This will prevent overheating of the treatment area, and reduce the possibility of any adverse effect.
- Keep 3 electrodes (**#2**) in good contact with skin.
- If the energy level feels uncomfortable, reduce it to a level that is comfortable for you.
- Only change the energy level when the device is not in contact with skin.

5.3 What to expect

During a treatment session it is normal to experience or feel:

- A sensation of warmth
- Some redness or slight swelling.
- this should subside within a few hours, or up to one day.
- Sound of the RF current.

5.4 Possible side effects

If you follow the instructions fully and pay attention to all warnings and contraindications, complications and side effects are uncommon. However, every cosmetic treatment, including those designed for home use, involves some degree of risk. Therefore it is important that you understand and accept the risks and complications that may occur.

- Excessive edema (swelling) and erythema (redness) are rare reactions to the treatment that, in most cases, subside within a few days.
- In case these signs appear, immediately use a damp cool cloth to cool the area for at least 10 minutes and apply an anti-burn cream.
- In rare cases the user may experience dilation of the superficial blood capillaries.
- This may result in a temporary redness for a few days.
- Also, some of the capillaries may present red spots that will also reduce with time.
- If these signs appear, we recommend that you cool the area with chilled water immediately.

- i** Consult your physician if you experience any persistent symptoms including pain, skin redness, swelling or a burning sensation.

6. Get started!

3. Make sure that your skin is clean and dry, and that you have removed all jewelry.
4. Remove the device and the other components from the box.
5. Start the device: Press the control button (**#5**). The device starts at energy level 1 automatically, and 1 energy level indicator light will turn on.

- i** Check whether the treatment medium has expired.

6. Apply a thin layer of treatment medium to the area you want to use.

- i** Treatment medium (or any other water-based medium) facilitates smooth movement of the treatment surface over the skin and enables conductivity. To ensure safe and efficient treatment, always start by applying medium and apply more medium if needed.

7. Using energy level 1, place the device on the area you want to treat. Make sure that it is in good contact with the skin. You will feel a mild sensation of warmth.
8. Move the device in gentle massage-like motions on the area you wish to treat. Progress slowly and carefully. Ensure 3 electrodes (**#2**) are in good contact with on the skin you want to treat when moving the device.
9. With the device constantly in motion, treat each area for at least 3-5 minutes.

10. Increase the energy level gradually. Using the highest energy level will get the best results, but it's important that you set the energy level that is comfortable for the area you are treating.

i You will feel the skin temperature increase as you increase the energy level. It usually takes about 1 minute for the device to reach the required temperature. During the treatment your skin will get a little red. Don't worry; this is normal.

 If you feel unwell, stop using the device immediately.

11. Try to move the device slowly at first. If the treatment area is hot and uncomfortable, try to move the device faster.
12. Clean the treated area: remove the medium residue off the treatment surface (**#1**), and clean the skin with warm water. Dry the skin and apply moisturizing facial products.

7. Cleaning, maintenance & storage

- Clean the device after each treatment.

 Never immerse your device or any of its parts in water!

- Turn off the device and put it in a power-off state
- Gently wipe off the care medium residue off the treatment surface (**#1**) of the device with damp paper (such as toilet paper or wipe).
- Make sure you remove all medium residue to prevent the medium from solidifying or damaging the device.
- Use a clean damp cloth to clean the device and the treatment surface (**#1**).

- You may use cleaning fluids with up to 65% alcohol, but never use acetone or other abrasive liquids, because they will damage your device.
- Finally, wipe the device and the treatment surface (#1) dry with a clean dry cloth.
- After cleaning, we recommend you store the device in its original box.
- Store the device in a cool and dry place.

8. Troubleshooting

Do not attempt to open or repair your device. You may expose yourself to dangerous electrical components and may cause serious harm to yourself. It will also void the warranty.

Problem	Checks
The device does not start.	Make sure the battery is not low.
I don't feel any heat from the device even though I pressed the control button.	Heat will only be emitted if the treatment surface (#1) is in good contact with your skin, and you have applied medium before you start your treatment.
One or more lights blink continuously.	The device has malfunctioned. Restart the device.
How to check product authenticity?	Visit www.silkn.hk and enter the product serial number.

For problems or queries, please contact Customer Service.

9. Customer service

For more information, please visit www.silkn.hk.

Contact Customer Service to report unexpected operations or events, or any other problem with your device. Our contact details are:

Hotline: 852-81045000

Email: info@fitboxx.com

10. Warranty

Please refer to the separate Product Warranty booklet for complete information.

- i** The manufacturer accepts no liability for damage or injuries caused by improper or incorrect use.

11. Disposal

If your device reaches its End of Life stage, be sure to follow local regulations for the disposal of electrical items. **Refer to local laws and regulations.**

- Dispose of the packaging separated into single type materials. Dispose of cardboard and carton as waste paper and film via the recyclable material collection service.
- Do not dispose the used components with household waste.
- If your device reaches its End of Life stage, be sure to follow your country's laws for the disposal of electrical items. Dispose in accordance with your country's laws and help protect the environment.

- Do not dispose with household waste! The WEEE directive requires that you dispose of this electrical device (including all accessories and parts) at a designated recycling collection point for electrical devices. This prevents the release of harmful substances* into the environment. Make sure that the device is in an uncharged state before disposal.

* Products labeled with Cd = cadmium, Hg = mercury, Pb lead.

12. Signs & symbols

	A statement or event that indicates a possible or impending danger, or problem when using the device.
	Important information about using or maintaining the device.
	Products marked with this symbol fulfil the requirements of the EU directives.
	Warning: Do not use in the bathroom.
	Do not dispose with household waste. Dispose in accordance with your country's laws and help protect the environment. In accordance with the WEEE Directive.
	The built-in rechargeable battery contains substances that may pollute the environment. Always remove the battery before you discard at an official collection point for batteries.

- This device is not suitable for use in the presence of flammable anaesthetic mixture with air or with oxygen or nitrous oxide.

13. Specifications

Technology	RF technology
LLLT	630nm±5nm
RF	1 MHz ±30%; Maximum 13W
Operation & safety	Heat stabilizer preventing skin overheating
System weight	50.9g
Temperature in operation	10°C to 35°C
Temperature in storage	-20°C to +60°C
Relative humidity in operation	30~75%
Relative humidity in storage	10~75%
Atmospheric pressure in storage	700 to 1060hPa
Atmospheric pressure in operation	500 to 1060hPa
Input power	DC 5V/1A
Battery capacity	500mAh
Charging method	Type - C charging port (#7)
Charging time	About 3 hours

14. Legal notice

Copyright 2021 © Home Skinovations Ltd. All rights reserved.

Home Skinovations Ltd reserves the right to make changes to its products or specifications to improve performance, reliability, or manufacturability. Information furnished by Home Skinovations Ltd is believed to be accurate and reliable at the time of publication. However, Home Skinovations Ltd assumes no responsibility for its use. No license is granted by its implication or otherwise under any patent or patent rights of Home Skinovations Ltd. No part of this document may be reproduced or transmitted in any form or by any means, electronic or mechanical, for any purpose, without the express written permission of Home Skinovations Ltd. Data is subject to change without notification. Home Skinovations Ltd has patents and pending patent applications, trademarks, copyrights, or other intellectual property rights covering subject matter in this document. The furnishing of this document does not give you any license to these patents, trademarks, copyrights, or other intellectual property rights except as expressly provided in any written agreement from Home Skinovations Ltd. Specifications are subject to change without notice. Silk'n and the Silk'n logo are registered trademarks of Home Skinovations Ltd, Tavor Building, Shaar Yokneam P.O. Box 533, Yokneam 2069206, ISRAEL, www.silkn.com, info@silkn.com.

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Copyright 2021 © Home Skinovations Ltd. 保留所有版权权利。